


Gas supply

- attach regulator to “zero air” bottle
- attach $\frac{1}{4}$ ” tubing to regulator (fasten with $\frac{9}{16}$ ” wrench)

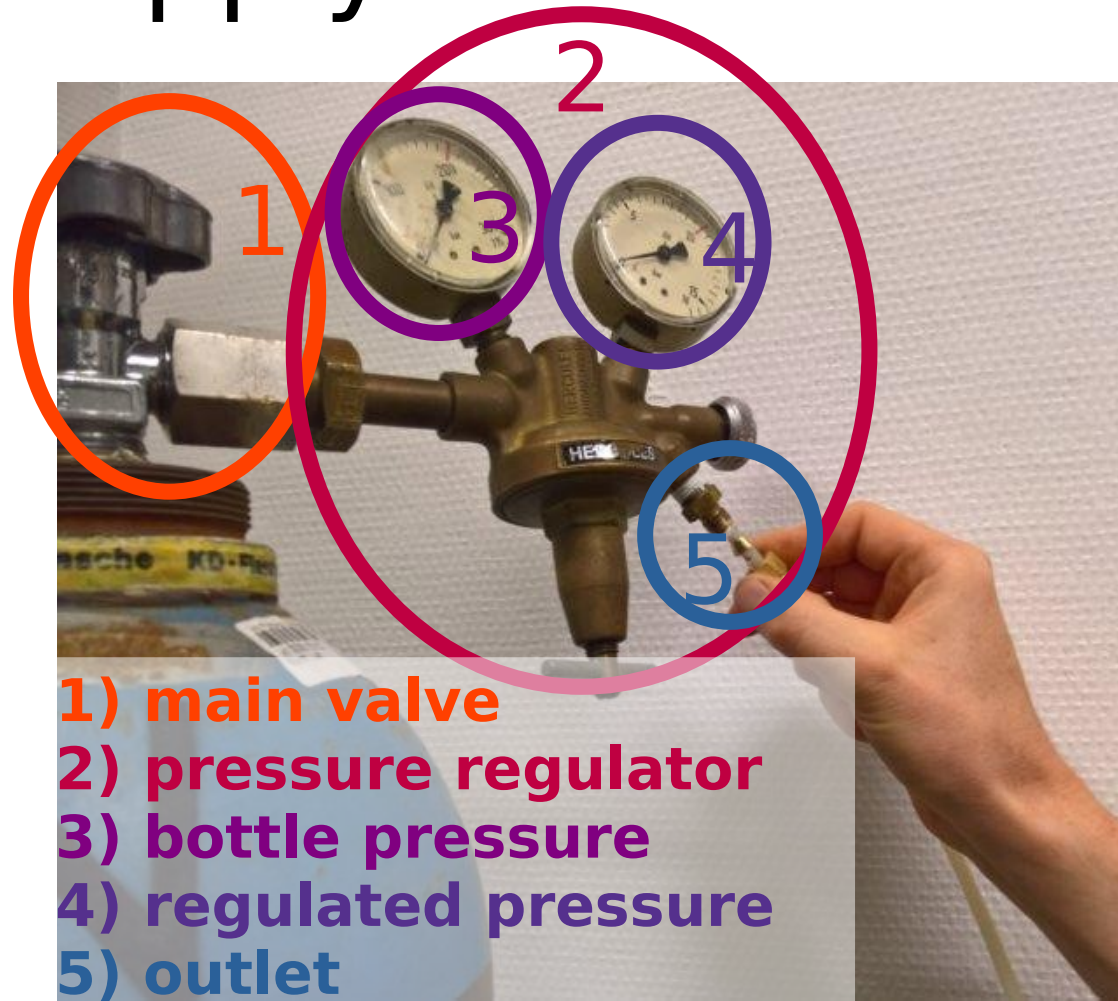


Gas supply

- The valve to regulate the pressure at the outlet might look different at different regulators



Gas supply

- use finger to block open side of $\frac{1}{4}$ " tubing
- open main valve
- regulate pressure (indicated by 4) to 1.5 to 2 bars
- close main valve

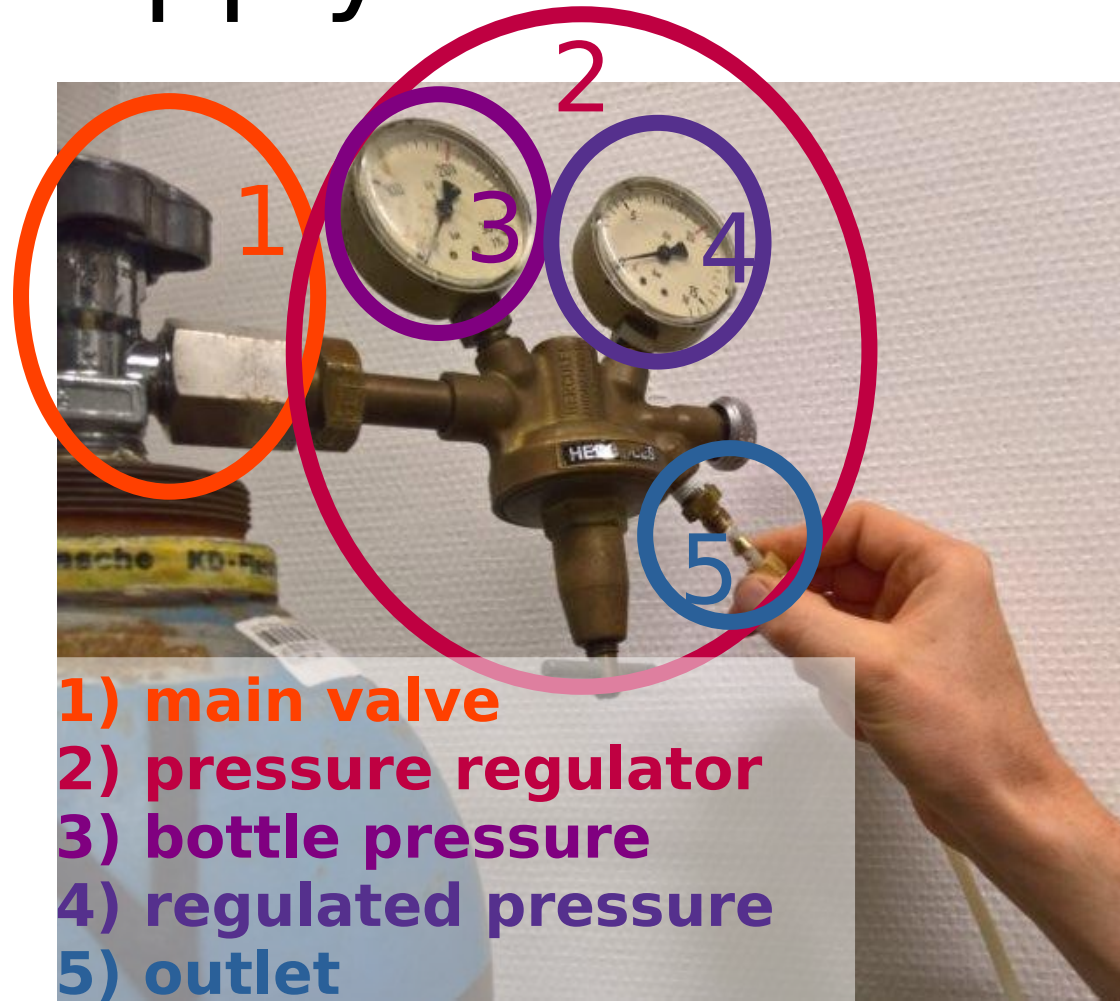


Mass flow controller inlets

- attach 1/4" tubing to T-connector
- attach T-connector to flow controller inlets




- fasten all connections with 9/16" wrench

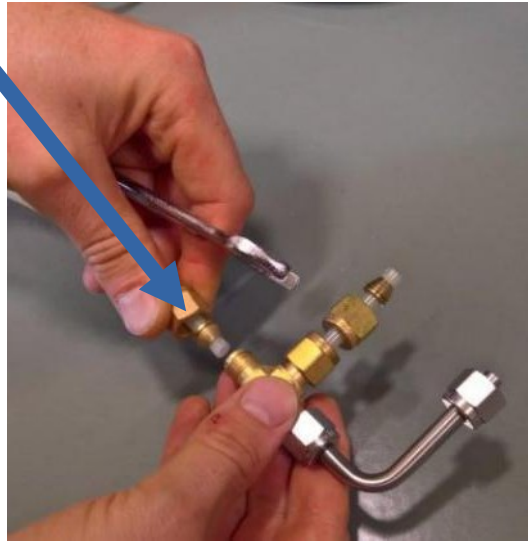

Mass flow controller outlets

- attach 1/8" tubing to flow controller outlets (black to black, red to red)
- fasten with 9/16" wrench

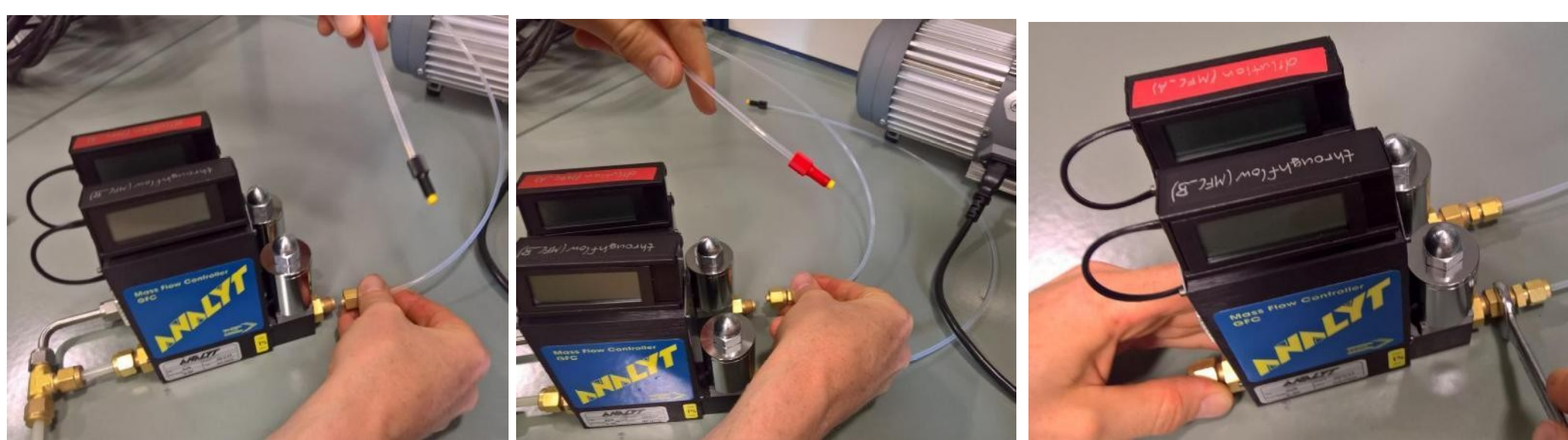




Picarro inlet

- attach 1/8" tubing terminated by blue screw to Picarro inlet
- fasten with 9/16" wrench



Plotbox manifolds

- screw the three 1/8" tubings whose other ends are 9/16" nuts (black & red @ flow controllers, blue @ Picarro) into sides of accordingly colored plotbox manifolds
- in case nut and ferrule are already pressed together, separate them before screwing

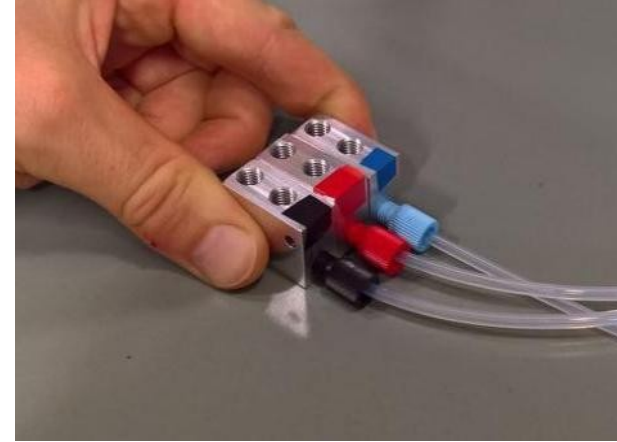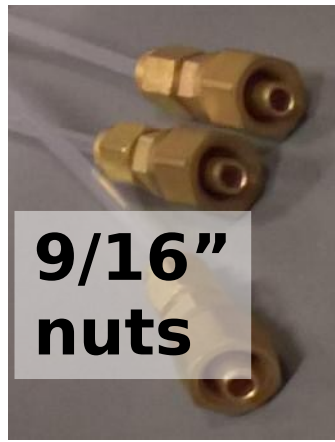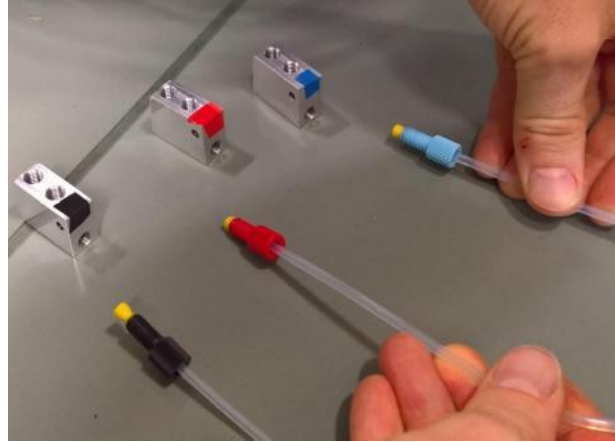


Plotboxes

- screw ends of plotbox connection tubing (3 x 1/8" with color coded screws on both ends) into plotbox manifolds and side ports of the valveblock manifolds of each plot box




- Tip: take valve unit out of plotbox for easier access

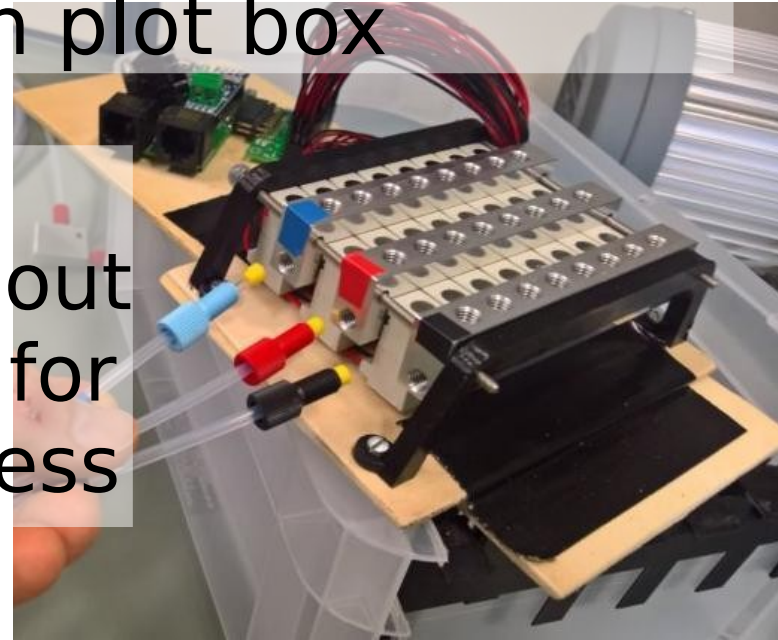


Plotboxes

- put valve unit back into the box


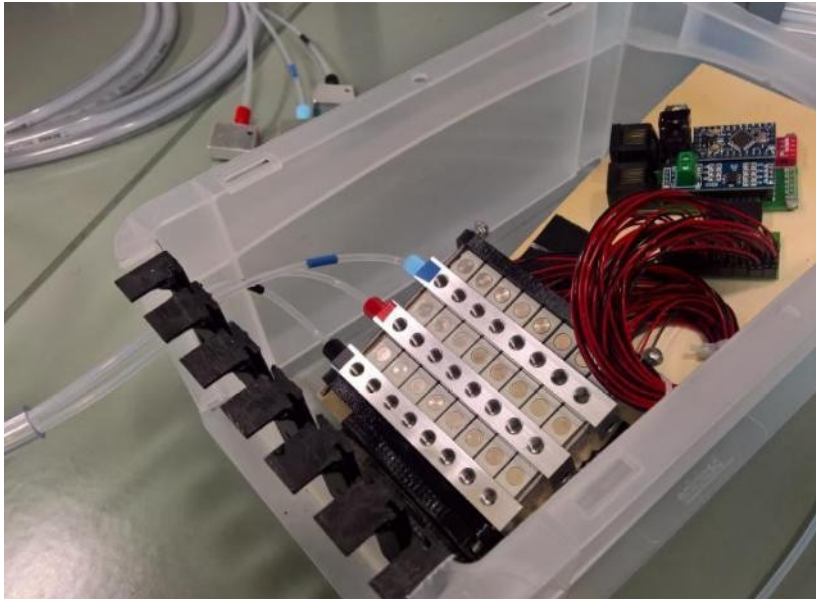


- guide tubing thorough leftmost access slot

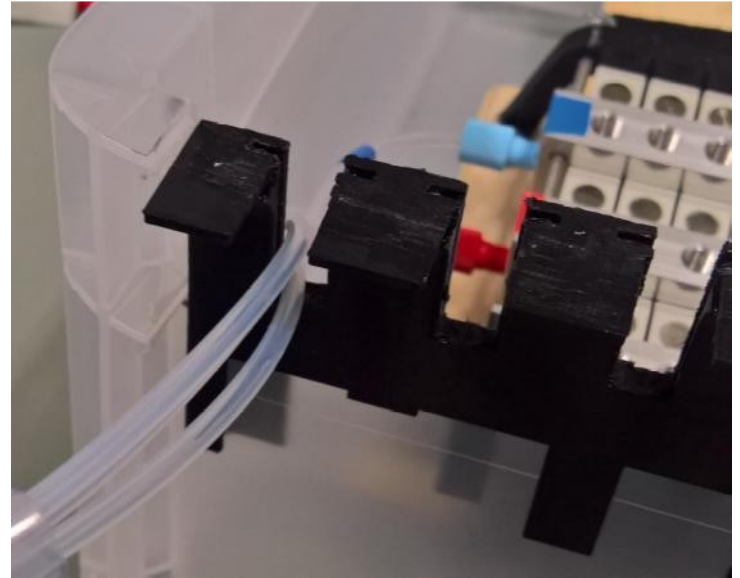


Probes

- make sure nut and ferrule are not pressed together
- screw all three lines of one probe into the same column (stick to the color coding)

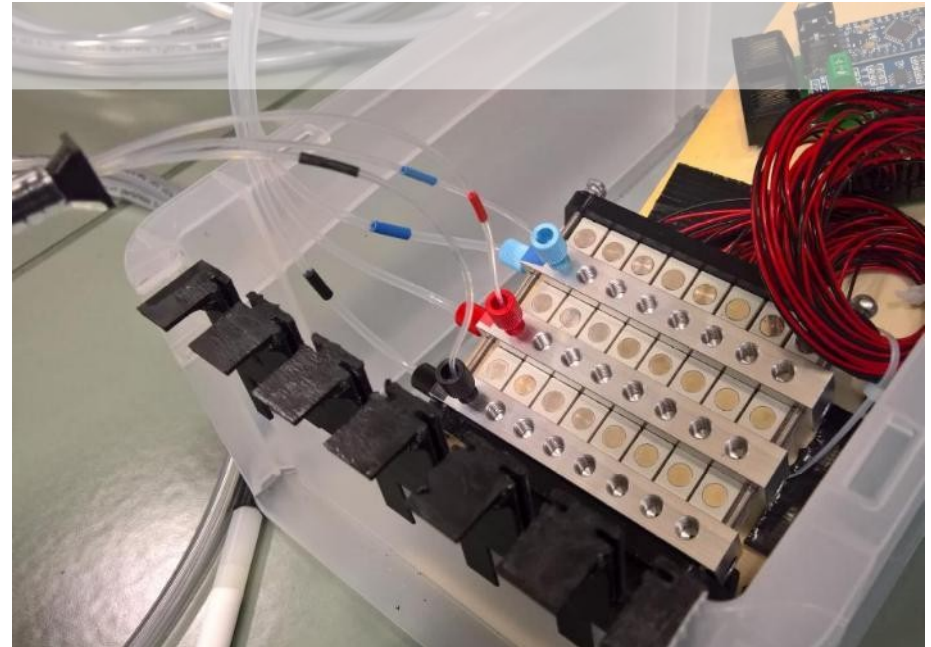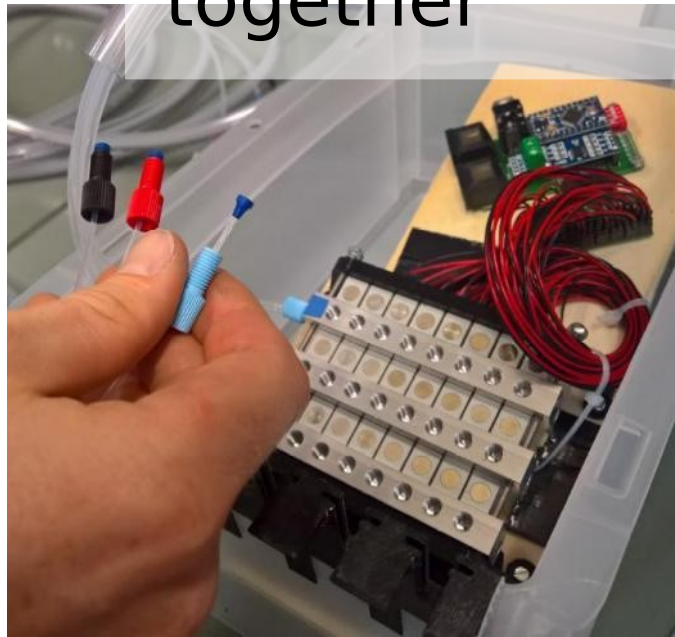

Probes

- slide end of protection tube into free slot of plotbox (one slot can take two probes)
- preferably use the four rightmost slots

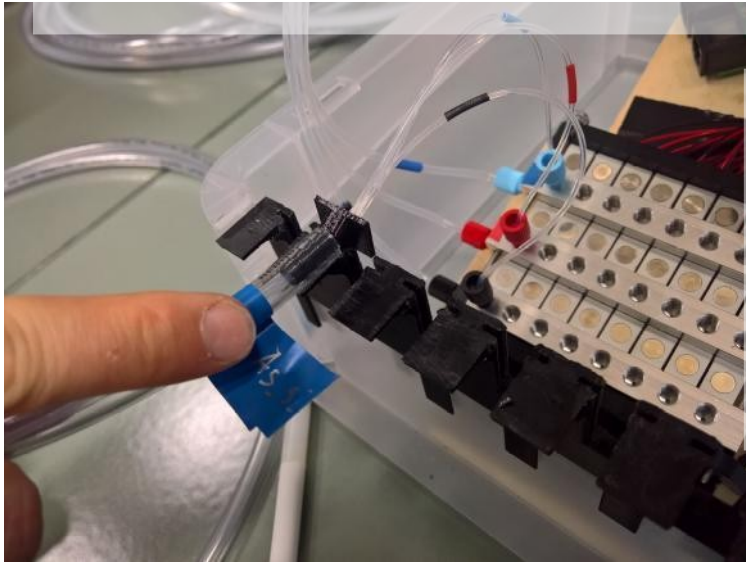

- repeat for remaining probes

